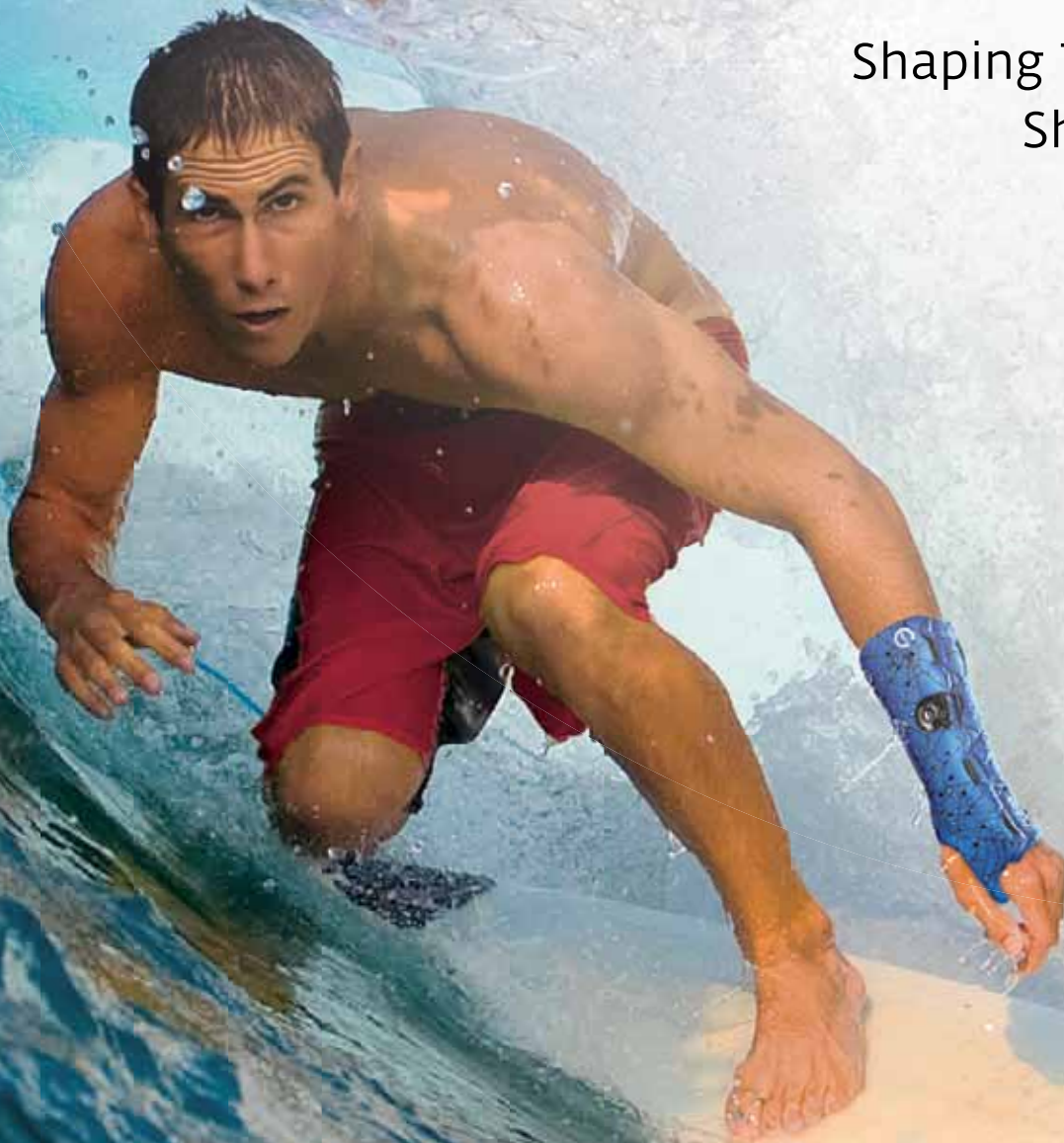




# Exos Bracing

Shaping Technology,  
Shaping Lives



# A Great Alternative to Conventional Casting!

Exos is the **only** bracing system offering you and your patients a removable, adjustable, reformable and waterproof solution for the treatment of fractures and other injuries requiring stabilization.

## Open it, Close it!

### REMOVABLE

Although the Exos brace is custom fit to your patient, it can be easily removed, allowing:

- Examination of the fracture site
- Examination of a wound or an incision
- Assessment of muscle atrophy
- The ability to initiate rehabilitation therapy earlier in the healing process
- Skin and brace cleaning
- Sawless removal

Comfortable, lightweight material

Exclusive BOA® Dynamic closure system

Antimicrobial treated

Washable and waterproof

BOA Locking Ring



# Warm it, Fit it!

## REFORMABLE

The Exos brace is completely reformable, allowing:

- Modification to accommodate anatomical changes as healing progresses
- Serial changes of position
- One brace throughout duration of treatment
- Less material waste

## ADJUSTABLE

Exos uses the exclusive BOA® reel and lace dynamic closure system, allowing:

- Expansion to allow for initial swelling and to reduce the risk of Compartment Syndrome
- Tension modification for custom fit despite atrophy or changes in soft tissue
- Improved patient comfort and fit



1. Simply heat brace in oven for 3-5 minutes



2. Open and slide onto patient's limb



3. Gently mold and adjust for a comfortable fit



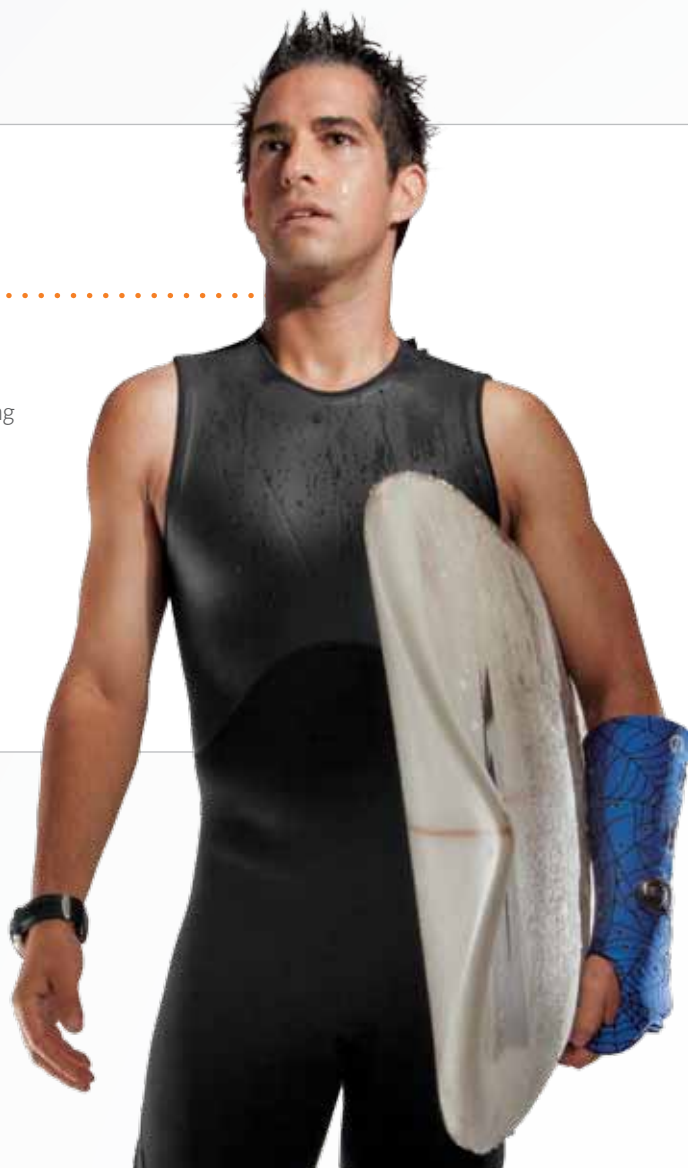
4. Use the reel and lace to tighten as needed

## It Takes to Water

---

The Exos brace is waterproof\*, allowing:

- Patients to continue their regular daily activities such as showering
- Children to remain active
- Athletes to continue water activities
- Use of blow dryer to quickly dry brace and skin



\*Limitations and precautions apply. See product instructions for use



## It's Easy to Clean

---

The Exos brace is designed for improved hygiene, allowing:

- Patients to hand wash and dry as needed
- Antimicrobial treated material to reduce odor
- Brace to be removed in order to clean skin or injury site



## It's See Through & Lightweight

The Exos brace is radiolucent, allowing:

- X-rays to be taken without the removal of the brace
- Clinical assessment of fracture healing without the need for re-casting
- The Exos braces are significantly lighter than plaster casts or fiber glass ones



## It's Original!

Exos braces come in a selection of different colors and patterns:

- Making them attractive and fun to wear, especially for children



1111 = Black



4444 = Pink



2293 = Spider



2388 = Camo



3285 = Polka dot



2265 = Denim

# Reformable, Removable Bracing System



### Wrist Brace with BOA

- Carpal bone & ligament injuries
- Stable wrist fractures



### Radial Gutter Fracture Brace

- 2nd & 3rd Metacarpal fractures



### Short Arm Fracture Brace

- Distal Radius Fractures (non-displaced)
- Wrist sprains, ligament injuries



### Long Thumb Spica with BOA

- Scaphoid fracture
- Radial or Ulnar wrist strain



### Short Arm Fracture Brace - Open Thumb

- Distal Radius Fractures (non-displaced)
- Wrist sprains, ligament injuries



### Boxer's Fracture Brace

- 4th & 5th Metacarpal fractures



### Thumb Spica Fracture Brace

- Skier's Thumb, Bennett's & Rolando's fractures
- Scaphoid fractures



### Short Thumb Spica

- Skier's Thumb, Sprains, Dislocations

### Extended Short Thumb Spica

- Thumb IP joint hyperextension, Pin fixation



### Hand based Ulnar Gutter Brace

- 4th & 5th Metacarpals or Phalanges fractures

### Hand Based Radial Gutter Brace

- 2nd & 3rd Metacarpals or Phalanges fractures

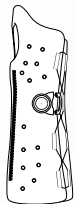


### Humeral Fracture Brace

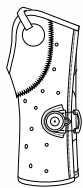
- Mid-shaft, non-displaced fractures

## ORDERING INFORMATION

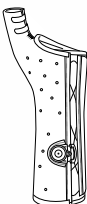
### SHORT ARM FRACTURE BRACE

Color options	Left	Right	Size
	N/A	N/A	XXS
	310-31	310-32	XS
	310-41	310-42	S
	310-51	310-52	M
	310-61	310-62	L
	310-71	310-72	XL

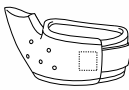
### WRIST BRACE WITH BOA

Color options	Left	Right	Size
	221-21	221-22	XXS
	221-31	221-32	XS
	221-41	221-42	S
	221-51	221-52	M
	221-61	221-62	L
	221-71	221-72	XL

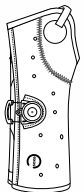
### THUMB SPICA FRACTURE BRACE

Black only	Left	Right	Size
	N/A	N/A	XXS
	315-31	315-32	XS
	315-41	315-42	S
	315-51	315-52	M
	315-61	315-62	L
	315-71	315-72	XL

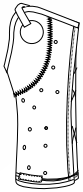
### SHORT THUMB SPICA

Color options	Left	Right	Size
	210-21	210-22	XXS
	210-31	210-32	XS
	210-41	210-42	S
	210-51	210-52	M
	210-61	210-62	L
	210-71	210-72	XL

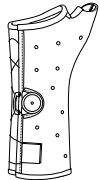
### SHORT ARM FRACTURE BRACE - OPEN THUMB

Color options	Left	Right	Size
	N/A	N/A	XXS
	312-31	312-32	XS
	312-41	312-42	S
	312-51	312-52	M
	312-61	312-62	L
	312-71	312-72	XL

### WRIST BRACE

Black only	Left	Right	Size
	N/A	N/A	XXS
	220-31	220-32	XS
	220-41	220-42	S
	220-51	220-52	M
	220-61	220-62	L
	220-72	220-72	XL

### LONG THUMB SPICA WITH BOA

Color options	Left	Right	Size
	231-21	231-22	XXS
	231-31	231-32	XS
	231-41	231-42	S
	231-51	231-52	M
	231-61	231-62	L
	231-71	231-72	XL

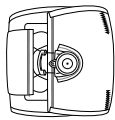
### EXTENDED SHORT THUMB SPICA

Black only	Left	Right	Size
	N/A	N/A	XXS
	211-31	211-32	XS
	211-41	211-42	S
	211-51	211-52	M
	211-61	211-62	L
	211-71	211-72	XL

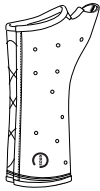
### PEDIATRIC SHORT ARM FRACTURE BRACE

Color options	Left	Right	Size
	311-11	311-12	XXS
	311-21	311-22	XS

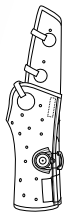
### HUMERAL FRACTURE BRACE

Black only	Left	Size
	330-40	Small
	330-50	Medium
	330-60	Large


### LONG THUMB SPICA

Black only	Left	Right	Size
	N/A	N/A	XXS
	230-31	230-32	XS
	230-41	230-42	S
	230-51	230-52	M
	230-61	230-62	L
	230-71	230-72	XL


### BOXER'S FRACTURE BRACE

Black only	Left	Right	Size
	N/A	N/A	XXS
	325-31	325-32	XS
	325-41	325-42	S
	325-51	325-52	M
	325-61	325-62	L
	325-71	325-72	XL


### HAND BASED ULNAR GUTTER BRACE

Black only	Left	Right	Size
	N/A	N/A	XXS
	125-31	125-32	XS
	125-41	125-42	S
	125-51	125-52	M
	125-61	125-62	L
	125-71	125-72	XL

### RADIAL GUTTER FRACTURE BRACE

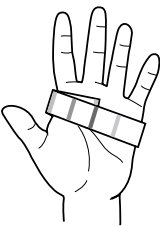
Black only	Left	Right	Size
	N/A	N/A	XXS
	326-31	326-32	XS
	326-41	326-42	S
	326-51	326-52	M
	326-61	326-62	L
	326-71	326-72	XL

### HAND BASED RADIAL GUTTER BRACE

Black only	Left	Right	Size
	N/A	N/A	XXS
	126-31	126-32	XS
	126-41	126-42	S
	126-51	126-52	M
	126-61	126-62	L
	126-71	126-72	XL

### SIZING CHART

Wrap Exos measuring tape around palm to measure size



Size	Centimeters	Inches
XXXS	14.5 - 16.5 cm	5¾" - <6½"
XXS	16.5 - 18.5 cm	6½" - <7¼"
XS	18.5 - 20 cm	7¼" - <8"
S	20 - 23 cm	8" - <9"
M	23 - 25.5 cm	9" - <10"
L	25.5 - 28 cm	10" - <11"
XL	> 28 cm	> 11"

DJO Global | 1a Guildford Business Park | Guildford | Surrey | GU2 8XG | ENGLAND

[www.DJOglobal.eu](http://www.DJOglobal.eu)



*Together in Motion™*